ESCO LIFESCIENCES





*Esco Powdermax™ 1 Powder Weighing Balance Enclosure Model PW1-3A_*

Powdermax™

Powder Weighing Balance Enclosure

Introduction

Laboratory operations involving weighing, mixing and dispensing hazardous powders can now be performed safely and economically with the Esco Powdermax™ 1 Powder Weighing Balance Enclosure. The enclosure's unique design maintains balance stability to 4 decimal places without compromising containment.

Designed specifically for powder containment applications, the Powdermax Enclosure employs carefully directed inward airflow at the face of the enclosure to protect the operator from exposure to hazardous powders or potent pharmaceutical compounds. The aerodynamically designed enclosure reduces airflow turbulence for maximum containment while minimizing the effect that air drafts may have on sensitive weighing balances. The enclosure design exceeds the containment requirements of ANSI/ASHRAE 110 as proven by independent type testing. All exhaust air is HEPA-filtered, removing virtually all particulate matter, before recirculation to the laboratory.

Designed for Enhanced Usability and Efficiency

The Powdermax™ is designed specifically to provide the operator with a high level of usability, comfort and visibility.

- A 13-degree sloped front allows easy access to the work zone, eliminating operator fatigue and increasing productivity.
- Transparent 6 mm (0.2") frameless acrylic front and sides enhance visibility and operator comfort during longer periods of operation. The self-supporting pre-tensioned hinges on the front window provide easy access during loading and start-up.
- Ergonomically designed apertures for hands simultaneously provide maximum movement within the work zone while protecting the operator. The curved front edge minimizes airflow turbulence and improves user comfort.
- Electronic ballast for the fluorescent lighting provides zero-flicker with increased energy efficiency, reliability and service life with a lower heat output.
- Two electrical pass throughs on the back wall of the unit provide convenient access to power sources for equipment.
- The spacious interior dimensions will easily accommodate large analytical balances.

Enhanced Filtration System


The enclosure's unique filtration system ensures the highest level of protection for the operator.

- The Powdermax™ provides enhanced safety to laboratory personnel by drawing air across the front opening in the hood at a controlled rate, preventing the escape of powders and particulate matters into the general laboratory environment.
- A unique filter mounting system eliminates bypass leaks inherent in conventional clamping systems.
- The Powdermax™ is equipped with a generously sized HEPA filter system, rated at >99.99% efficiency for 0.3 micron sized particles.
- The aerodynamically designed system constructed of 1.2 mm (18 gauge) epoxy powder-coated electrogalvanized steel delivers maximum containment by ensuring airflow uniformity throughout the work chamber of the cabinet.
- Since HEPA-filtered air is recirculated to the laboratory, the enclosure does not have to be connected to an extraction/ exhaust system, resulting in low operating costs and mobility.

User-Friendly Control System

The user-friendly Esco Sentinel™ Silver microprocessor based control system is fully configurable according to operator's requirements and comes equipped with a number of enhanced features to promote cabinet usability.

Filtration System

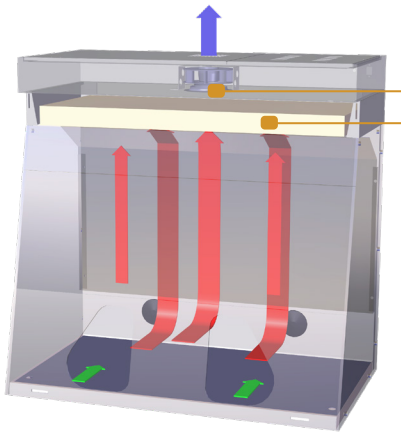

- Room air is drawn in across the front of the cabinet with an average velocity of 0.5 m/s (100 fpm).
- Air is drawn up through the cabinet's work zone and forced through the HEPA filter (>99.99% typical efficiency for 0.3 micron sized particles).
- The state-of-the-art baffle system ensures airflow uniformity throughout the cabinet's main chamber.
- The HEPA filtered air then returns to the laboratory stripped of all airborne contaminants or is vented through the optional exhaust collar to exhaust ducting for enhanced safety.

HEPA-filtered air Unfiltered / potentially contaminated air Inflow air

- Accurate true airflow velocity sensing technology measures all critical cabinet airflow parameters allowing superior monitoring. Temperature compensated sensors ensure increased accuracy.
- Password-protected administration can be set to restrict access to the main menu.
- A bright, easy-to-read, LCD display provides continuous monitoring of cabinet airflow.
- Audible and visual alarms ensure operator and environmental protection by alerting the user in the event of low airflow or unsafe sash positions.
- When the cabinet is switched on, its automatic warm-up cycle is activated, purging all contaminants from the work zone prior to use.
- An automatic post-purge cycle can be configured on shutdown, purging all residual contaminants before deactivation.

Highest Quality Enclosure Construction

All Esco cabinets are constructed to the highest quality using the finest materials.

- The enclosure's head unit is constructed of electrogalvanized steel with abrasion-resistant powder-coated finish.
- External surfaces are coated with Esco Isocide™ antimicrobial coating to protect against surface contamination and inhibit bacterial growth. Isocide™ eliminates 99.9% of surface bacteria within 24 hours of exposure.
- The enclosure is supplied with a chemical and abrasion resistant phenolic resin work surface. The solid work surface minimizes transmission of fan vibration to the user's weighing balance.

Blower Efficiency

The Powdermax™ incorporates an energy efficient external rotor motor design which reduces operating costs and has extremely low noise and vibration levels. The cabinet blower system creates an inflow of air from the ambient environment into the cabinet with an average velocity of 0.5 m/s (100 fpm).

Options and Accessories

The Powdermax™ is available with a number of options and accessories to meet your needs. These include:


- Mobile cart with foldable trays
- Height-adjustable lab chair
- Ergonomic foot rest


Further customization specifications are available upon request.

Microbalance Stability Testing

Weighing Accuracy Test on Esco Powdermax™ Powder Weighing Balance Enclosure

The Esco PowderMax™ cabinet at nominal airflow can provide a stable work condition for sensitive powder weighing work with extremely low standard deviation in the range of 0.000001 to 0.000011 for weighing 20g, 5g, and even 0.5g objects.

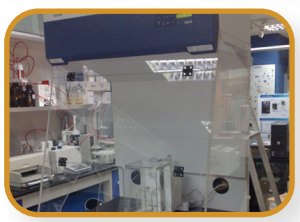

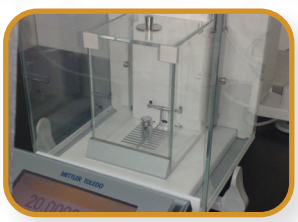

Surrogate Powder Testing


Active Pharmaceutical Ingredients (API) Surrogate Containment Testing on Esco Powdermax™ Powder Weighing Balance Enclosure


The Esco PowderMax™ Powder Weighing Balance Enclosure, operating at 0.50 m/s (100fpm) average inflow, was demonstrated to contain surrogate API (naproxen sodium), at an OEL of 26.3 nanogram / m^3, for a challenge of repetitive weighing of 20 x 1-gram quantities of API, independent of the operator's level of skill.



General Specifications, Powder Weighing Balance Enclosure							
Model		**PW1-3A1 2160001**	**PW1-3A2 2160003**	**PW1-3A3 2160004**	**PW1-4A1 2160010**	**PW1-4A2 2160006**	**PW1-4A3 2160007**
Nominal Size		0.9 meters (3 ft)			1.2 meters (4 ft)		
External Dimensions (W x D x H)		840 mm x 700 mm x 1125 mm (33.0" x 27.6" x 44.3")			1140 mm x 700 mm x 1125 mm (44.9" x 27.6" x 44.3")		
Internal Work Area, Dimensions (W x D x H)		820 mm x 650 mm x 829 mm (32.3" x 25.6" x 32.6")			1120 mm x 650 mm x 848 mm (44.1" x 25.6" x 33.4")		
Standard Filtration Elements	Main Filter	HEPA filter, typical efficiency of >99.99% at 0.3 microns, removes particles and aerosols					
Control System		Sentinel™ Silver Microprocessor Controller					
Inflow Air Velocity		Initial setpoint: average of 0.5 m/s (100 fpm)					
Sound Emission		<55 dBA at initial blower speed setting measured at typical operator work position					
Cabinet Construction	Main Body	1.2 mm (0.05") 18 gauge electrogalvanized steel with white oven-baked epoxy-polyester Isocide™ antimicrobial powder coated finish					
	Back Wall						
	Front Window	6 mm (0.2") Acrylic					
	Side Walls						
	Work Top	Phenolic Resin					
LED Lamp Intensity at Zero Ambient		>600 lux (>56 ft. candles)					
Electrical	Electrical Code	220-240 VAC, 50 Hz	110-120 VAC, 60 Hz	220-240 VAC, 60 Hz	220-240 VAC, 50 Hz	110-120 VAC, 60 Hz	220-240 VAC, 60 Hz
	FLA (A)	0.5	1	0.5	1	1	0.5
	Nominal Power (W)	78	87	85	87	81	87
Sash Specification	Sloped Front	13°					
	Sash Type	Hinged					
Net Weight		110 kg (242.5 lbs)			160 kg (352.7 lbs)		
Shipping Weight		120 kg (264.6 lbs)			170 kg (374.8 lbs)		
Shipping Dimensions, Maximum (W x D x H)		1050 mm x 1100 mm x 720 mm (2314.9" x 2425.1" x 1587.33")			1300 mm x 1100 mm x 720 mm (2866.0" x 2425.1" x 1587.33")		
Shipping Volume, Maximum		0.83 m^3 (29.3 cu.ft.)			1.0 m^3 (35.3 cu.ft.)		

Powdermax™ Powder Weighing Balance Enclosure Engineering Drawing

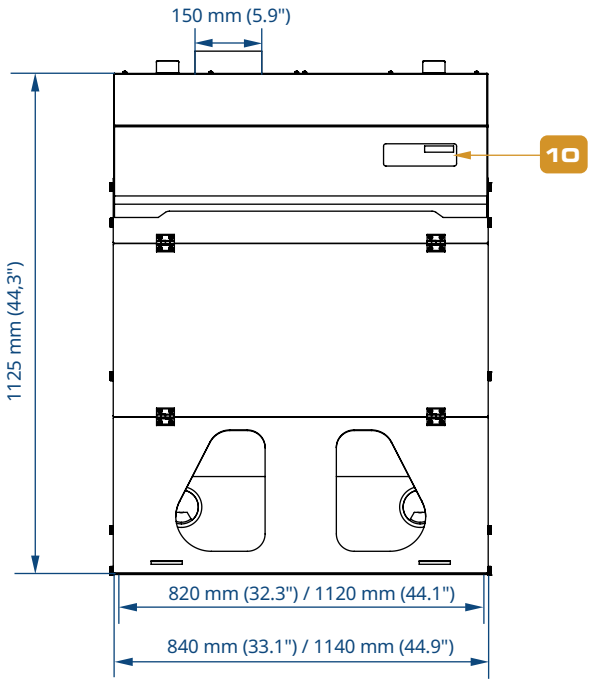

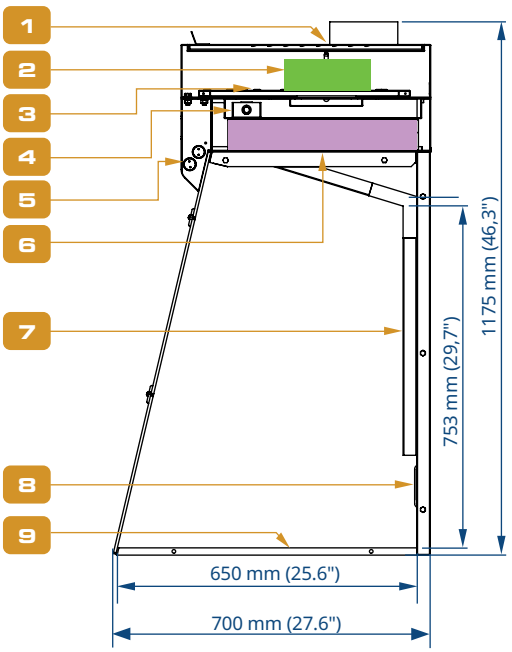


1. *Exhaust Collar*
2. *Fan Blower*
3. *Electrical Panel*
4. *Airflow Sensor*
5. *LED Lamp*
6. *HEPA Filter*
7. *Baffles*
8. *Pass-through ports with cover*
9. *Phenolic Work top*
10. *Esco Sentinel™/Silver Microprocessor control System*

Standards Compliance	Containment	Electrical Safety
	ANSI/ ASHRAE 110-1995, USA	UL-61010-1, USA CAN/CSA-C22-2 No.61010-1 EN-61010-1, Europe IEC61010-1, Worldwide

ESCO LIFESCIENCES GROUP

42 LOCATIONS IN 21 COUNTRIES ALL OVER THE WORLD

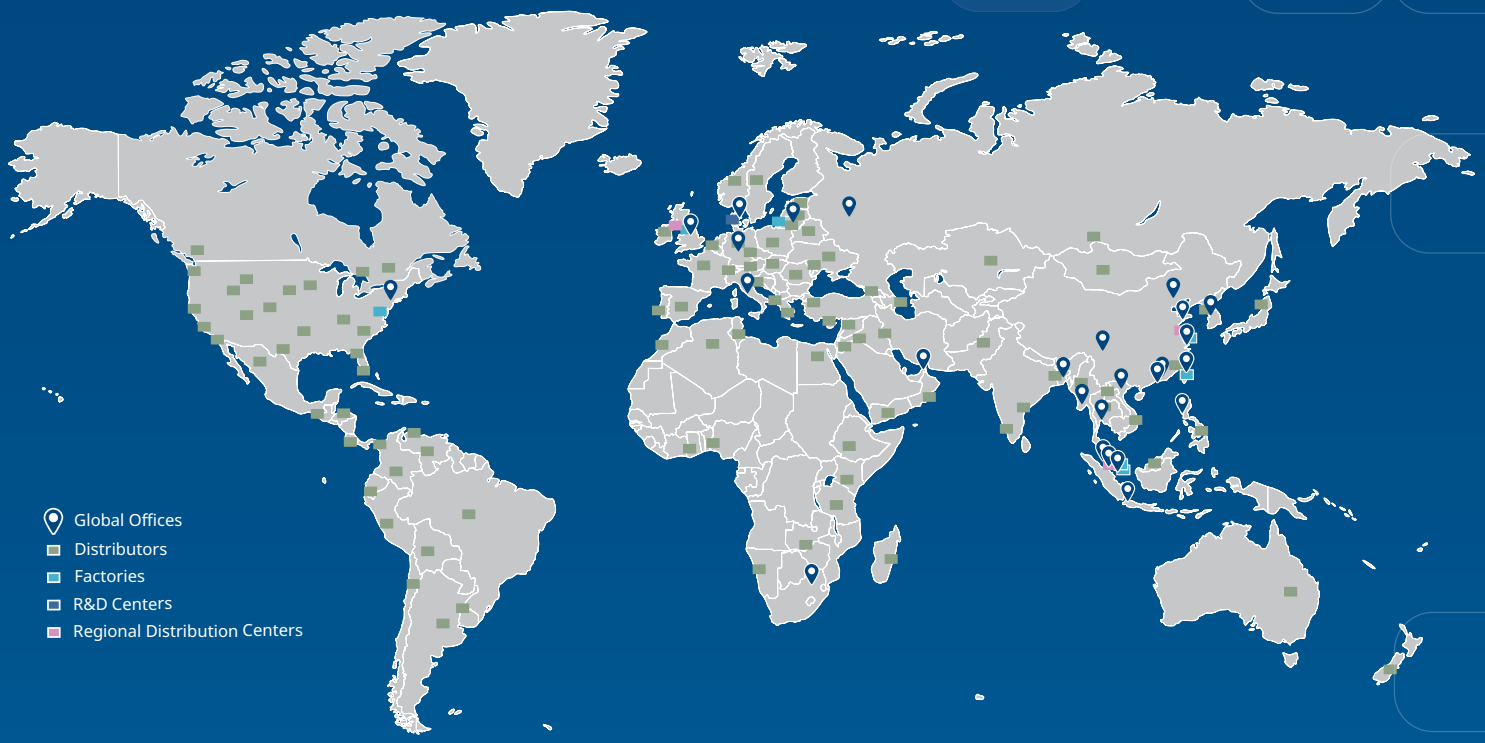

Follow us on social media, download our apps, and scan the QR code for more info.

@EscoLifesciences | @EscoLifesciences | @EscoLifesci | @Esco | @EscoLifesciences | @EscoLifesciences | Google Play Esco Lifesciences | Google Play Esco Lifesciences

ESCO® LIFESCIENCES GROUP

Esco Micro Pte. Ltd. • 19 Changi South Street 1 • Singapore 486 779
Tel +65 6542 0833 • mail@escolifesciences.com
www.escolifesciences.com

Esco Technologies, Inc. • 903 Sheehy Drive, Suite F, Horsham, PA 19044, USA
Tel: +1 215-441-9661 • Fax 484-698-7757
eti.admin@escolifesciences.com

Esco Lifesciences Group Offices: Bangladesh | China | Denmark | Germany | Hong Kong | India | Indonesia | Italy | Japan | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam

9010060_Powder Weighing_Balance_PW1-A_SellSheet_A4_vI_26
Esco can accept no responsibility for possible errors in catalogues, brochures and other printed materials. Esco reserves the right to alter its products and specifications without notice. All trademarks and logotypes in this material are the property of Esco and the respective companies.